

ELLIPSE™

Occipito-Cervico-Thoracic
Stabilization System


Comprehensive Implant Offering



Non-Threaded Locking Cap



ElliptiClick™ Rod Retention

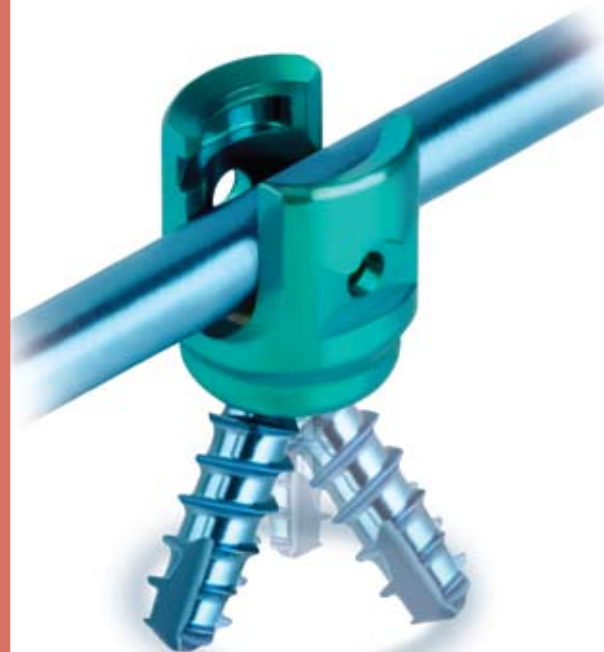

80° Total Angulation

ELLIPSE™ Stabilization System Benefits

Non-Threaded Locking Cap

- Eliminates cross-threading
- 90° rotation captures rod
- Set screw locks screw angle
- Low torque locking mechanism



Top-Loading T-Connector

- One set screw tightens the clamp onto the rod and secures the medial-lateral position
- Multiple implant sizes and sliding clamps eliminate the step of cutting the rod to size
- Clamp snaps over rod for controlled insertion
- Friction bushing holds the adjusted clamp position during insertion
- Adjusts to avoid contact with the dura or lower profile
- Straight and curved to accommodate patient anatomy



Head-to-Head T-Connector*

- Connects to polyaxial screw heads
- Allows insertion when screws are immediately adjacent
- Curved design accommodates patient anatomy
- Available in seven sizes

Parallel Connectors

- Allows attachment of ELLIPSE™ to titanium rod systems from 3.5mm to 6.5mm in diameter



Polyaxial Screw*

- 80° screw angulation ($\pm 40^\circ$) provides intra-operative flexibility
- ElliptiClick™ feature retains rod once inserted into screw head for easier cap insertion
- Elliptical shape provides additional resistance against screw head splaying

Lateral Connector

- Utilizes ElliptiClick™ feature for easier insertion
- May be cut to the appropriate length



Hooks

- Utilizes ElliptiClick™ feature for easier insertion
- Straight, left and right hook options
- 5mm and 7mm Offset Hooks



Rods

- 3.5mm diameter
- Tapered rod transitions from a diameter of 3.5mm to 4.0, 4.5, 5.0, 5.5, 6.0, 6.35, 6.5mm

Occipital Implants

- Large and small plates with seven points of fixation
- Pre-contoured plate with bending zones
- Rod acceptors allow for translation and rotation prior to final tightening

All ELLIPSE™ implants are made of titanium alloy

*For use in the thoracic spine (T1-T3) only



www.globusmedical.com

Globus Medical
Valley Forge Business Center
2560 General Armistead Avenue
Audubon, PA 19403

EC REP: RMS - UK Limited
28 Trinity Road, Nailsea, Somerset, BS484NU England

CE 0297

Customer Service:
Phone 1-866-GLOBUS1 (1-866-456-2871)
Fax 1-866-GLOBUS3 (1-866-456-2873)

©2009 Globus Medical. All rights reserved. ELLIPSE is a trademark of Globus Medical. Please refer to the product insert for description, indications, contraindications, warnings, precautions, and other important information.

GMSPS62
06.09