


# PROVIDENCE™ Anterior Cervical Plate System




Variable Angle Screws

Integrated Blocking Set Screw

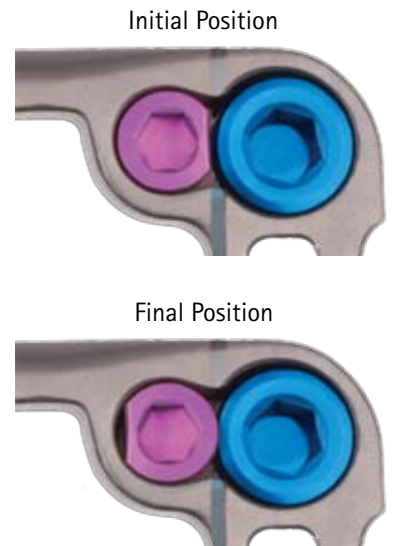

Audible, tactile and visual confirmation of screw blocking

DTS Guide, Pre-Set Angle Plate Holder

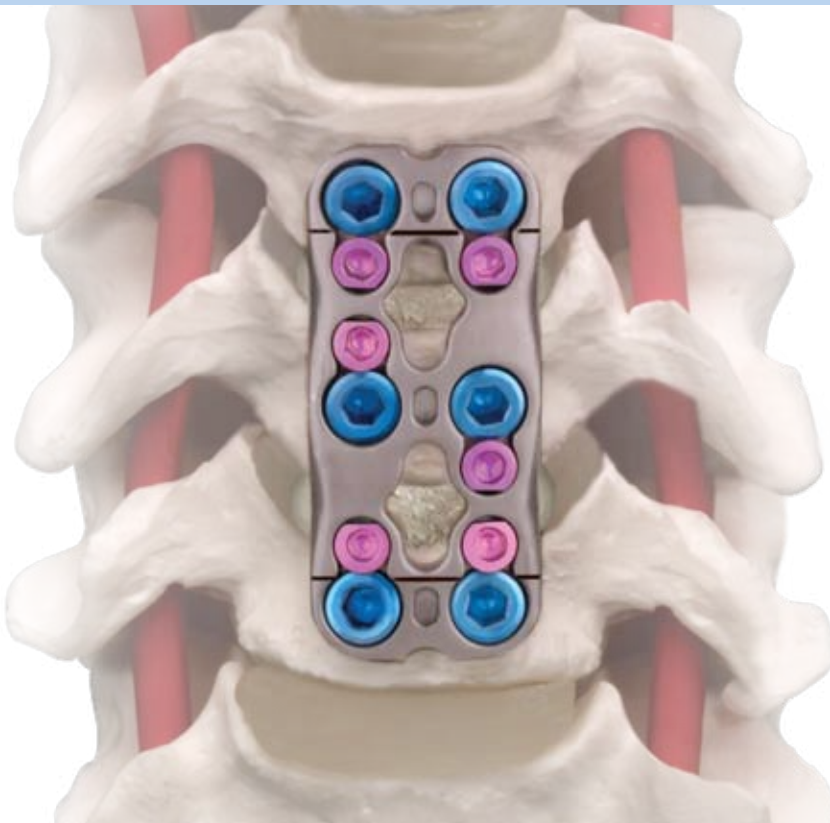

# PROVIDENCE™ Anterior Cervical Plate System Features

## Blocking Set Screw



- Provides audible, tactile and visual confirmation of blocking
- Pre-assembled into the plate
- Prevents screw backout

## Plates

- 2.3mm low-profile
- 16mm width
- Large graft visualization windows
- 1, 2, 3, 4 and 5-level plates
- Lengths from 10mm to 104mm


## Variety of Screw Sizes

- 4.2mm or 4.6mm diameter
- Self-Drilling or Self-Tapping
- Lengths from 10mm to 26mm
- Variable and Fixed Angle Screws

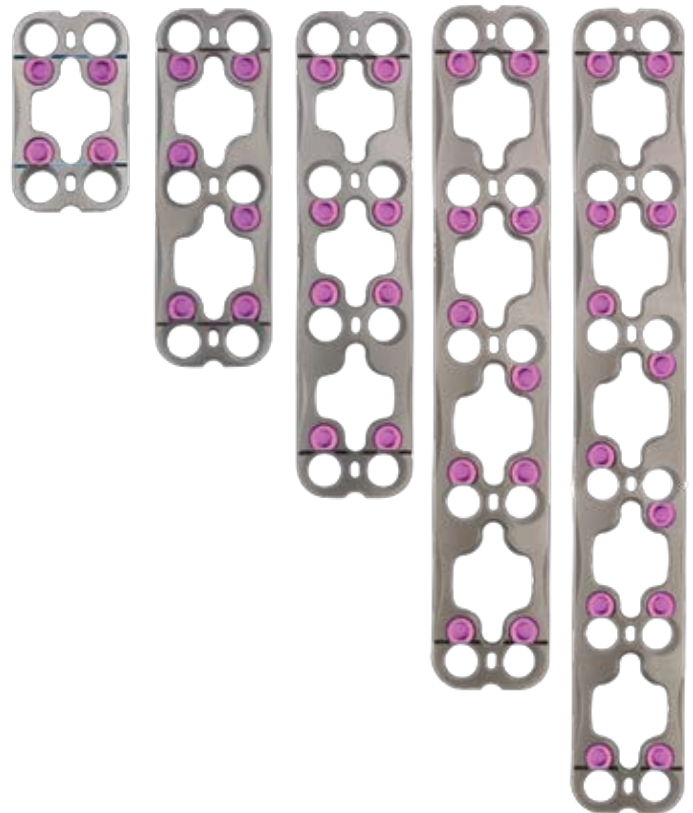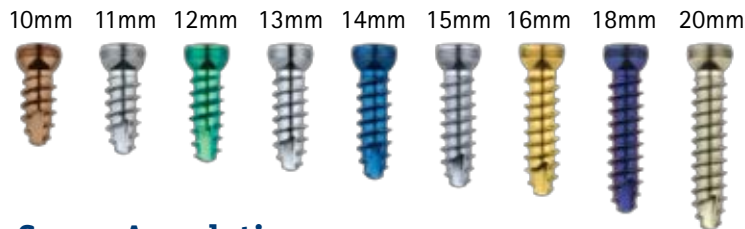


Variable Angle Screw

Fixed Angle Screw

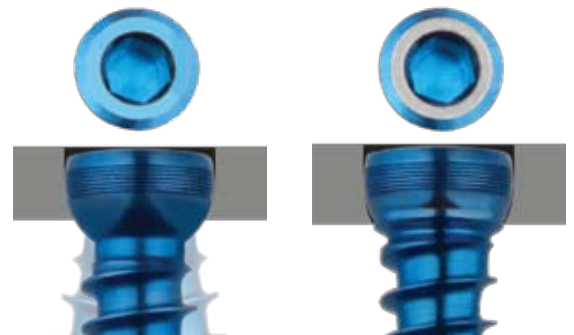

## Screw Angulation

- Variable Angle Screws provide +/- 20° angulation
- Pre-set angulation with drill guide
  - 6° medial
  - 10° cephalad/caudal respectively

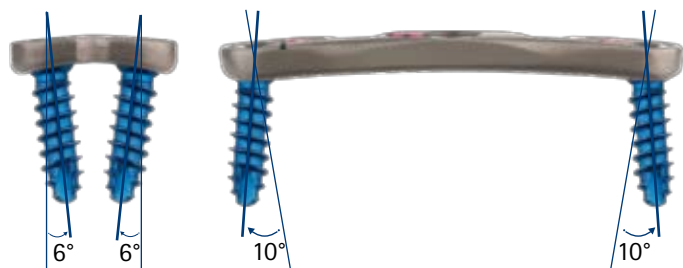

[www.globusmedical.com](http://www.globusmedical.com)

Globus Medical  
Valley Forge Business Center  
2560 General Armistead Avenue  
Audubon, PA 19403

Customer Service:  
Phone 1-866-GLOBUS1 (1-866-456-2871)  
Fax 1-866-GLOBUS3 (1-866-456-2873)

CE REP: RMS - UK Limited  
28 Trinity Road, Nailsea, Somerset, BS484NU England



©2009 Globus Medical. All rights reserved. PROVIDENCE is a trademark of Globus Medical. PRESERVE is a registered trademark of Globus Medical. Please refer to product insert for description, indications, contraindications, warnings, precautions, and other important information.